


H/S Elliptosphere™ Catheter Set and Tray

The Safe and Simple Answer for Both
Sonohysterography and Hysterosalpingography

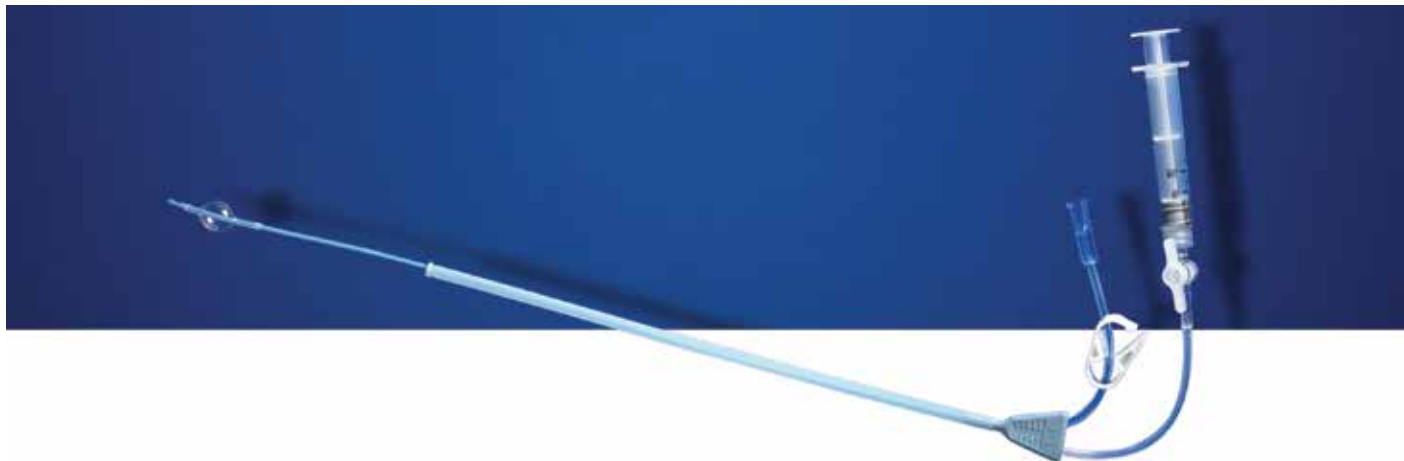

The patented Elliptosphere H/S catheter offers several advantages over the traditional balloon catheters including the ability to gradually inflate the balloon for placement in any size cervical canal, the ability to see the entire lower uterine segment during an SIS procedure, and greater patient comfort.

The Elliptosphere balloon contour and variable size permit secure placement in the cervical canal or uterine cavity for optimal visualization. The balloon's unique design prevents retrograde leakage of distension media, which facilitates peak uterine distention. A small catheter size enables placement next to ultrasound probe for minimal patient discomfort. The H/S Elliptosphere Catheter can be used with oil- or saline-based medium. Available in two versions – with and without an integrated stylet – the H/S Elliptosphere features latex-free construction.



Item #	Description
61-4005	H/S Elliptosphere Catheter, 5 Fr (10/Box)
61-4005S	H/S Elliptosphere Catheter Set with Integrated Stylet (10/Box)



DISTRIBUTED BY:
ADVANCED SURGI-PHARM INC.
www.advancedsurgipharma.com
Tel: 1-800-661-5432

Ver.04-22-2021