


Irrijet®

Splatter-Gard™

Simple, one-hand, thumb-controlled pressure provides an appropriately forceful stream (8 psi) directly to the desired area (or a gentler stream if needed). The soft flexible tip has an inside diameter equivalent to an 18-19 gauge needle.

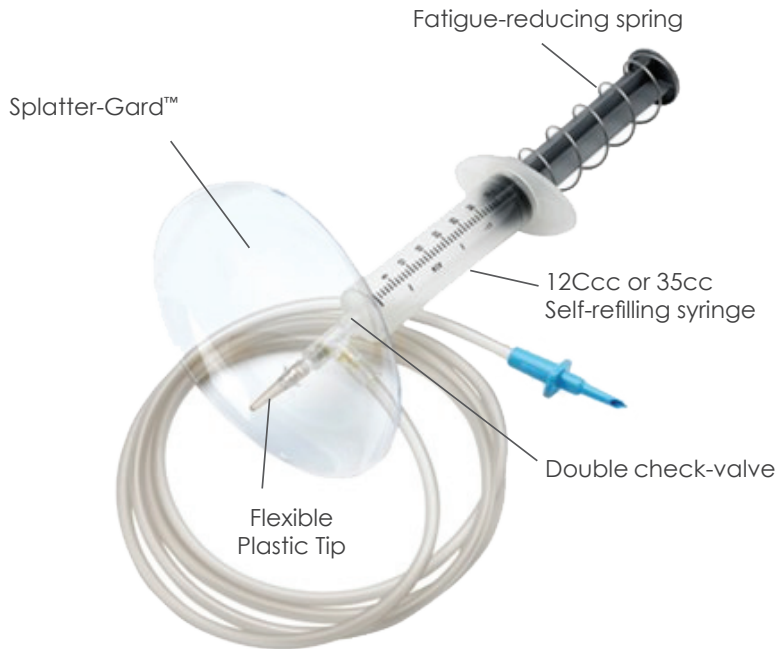


- Splatter-Gard™ prevents contamination from body secretions – new concave design provides maximum splash-back protection
- Soft, flexible plastic tip eliminates needle-stick for user, minimizes trauma to patient
- Irrigant is delivered with force equivalent up to 8 psi
- Latex-free
- Sterile

Item #	Description
55-0004	Irrijet GP (5/Box)
55-0504	Multi-Use Irrijet (5/Box)
55-2004	Single-Use Irrijet Z (5/Box)
55-2104	Single-Use Irrijet (5/Box)
55-2134	Single-Use Irrijet (5/Box)
55-6504	Irrijet DS (5/Box)
59-5014	Infuset GP (5/Box)

Irrijet GP

55-0004


Spring-loaded syringe connected to a double check-valve, through a side arm and 6-ft. length of tubing to a spike for entry into a bag of irrigating solution. Double check-valve prevents retrograde flow. Alternate movement of the syringe plunger delivers irrigating solution as needed. Syringe Capacity: 12cc.

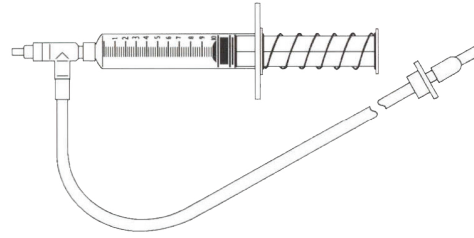

Multi-Use Irrijet

55-0504


Same basic assembly as 55-0004, plus one package of five individually sealed jet tips.

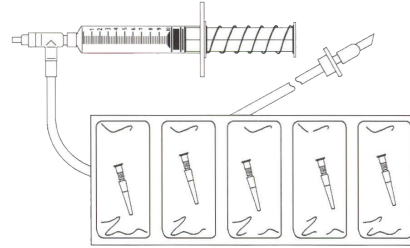

Single-Use Irrijet Z

55-2004

Same basic assembly as 55-0004, plus one Zerowet Splashshield®.

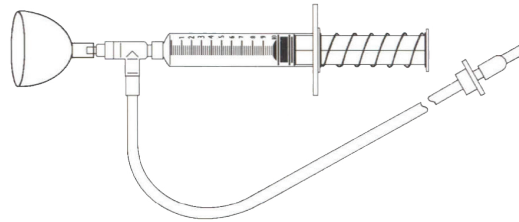

Single-Use Irrijet

55-2104

Same basic assembly as 55-0004, plus includes permanent flexible plastic jet tip and Splatter-Gard™.

55-2134

Same basic assembly as 55-2104, but with 35 cc capacity syringe.

