

SCOPE TRAYS



Scope Trays provide you with the protection you need for your expensive endoscopes and cameras.

Scope Trays are available in a variety of popular sizes to satisfy most of your instrument requirements. The trays are designed to withstand sterilization protocols for steam, ethylene oxide, dry heat, or cold sterilizing solutions. All instrument trays are manufactured with durable medical grade polymer.

Instrument Trays are ideal for storing, transporting, and protecting surgical medical instruments. Easy customization can be achieved for a variety of special configurations utilizing the unique Plug 'n Sert silicone bar system. This valuable feature allows you to place the silicone bar exactly where you need it for the most secure positioning.

Customization

Trays can be customized for individual needs.

Protection

Trays protect your instruments from costly damage.

Quality / Durability

Trays are manufactured with strict quality standards and are engineered to last.

Versatility


Trays can be used with a wide variety of instruments and can be individually customized.

Cost Savings

Unique Gridlock System™ allows for the use of one tray for many applications.

Sterility Assurance

The raised mat assures proper penetration of sterilizing agents between the mat and the tray surface.

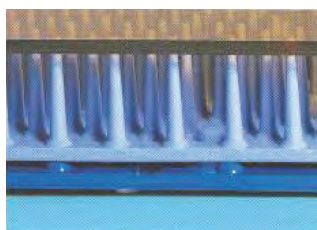

Trays are engineered with a unique "pebble like" interior surface. This raises the mat to assure proper penetration of the sterilization agent between the mat and the tray surface. Independent laboratory test have shown that the space between the mat and the surface is effectively penetrated.



Our latch is designed to last! Compare the quality and manufacturing standards of engineered trays and see the difference. Mold flow analysis gives products the strength and durability to satisfy your requirements.

Part no.	Description	Dimensions	
		Inches (WxLxD)	Cm (WxLxD)
SCOPE TRAYS 6500 SERIES			
TR-920-RO	Scope Tray - Base, Lid, & Bars	3.0 x 22 x 1.5	7.6 x 55.9 x 3.8
INSTRUMENT - MINI TRAYS 1300 SERIES			
TR-100-RO	Mini Tray - Base, Lid, & Bars	1.5 x 2.7 x 0.75	3.8 x 6.8 x 1.9
DIAMOND KNIFE TRAYS - 2500 SERIES			
TR-411-RO	Bar 2 Slots	-	-
TR-431-RO	Bar 3 Slits	-	-
TR-432-RO	Bar 6 Slits	-	-
INSTRUMENT - MINI TRAYS 2600 SERIES			
TR-500-RO	Micro Tray - Base, Lid & Mat	2.5 x 6.0 x 0.75	6.4 x 15.2 x 1.9
TR-510-RO	Micro Tray - Base, Lid & Mat	2.5 x 6.0 x 1.25	6.4 x 15.2 x 3.2
MICROSURGICAL INSTRUMENT TRAYS - 2700 SERIES			
TR-520-RO	Micro-Tray - Base, Lid & Mat	2.5 x 7.5 x 0.75	6.3 x 19 x 1.9
TR-521-RO	Micro-Tray - Base & Lid	2.5 x 7.5 x 0.75	6.3 x 19 x 1.9
MICROSURGICAL INSTRUMENT TRAYS - 4700 SERIES			
TR-701-RO	Micro Tray - Base & Lid	4.0 X 7.5 X 0.75	10.1 x 19 x 1.9
TR-710-RO	Micro Tray - Base, Lid & Mat	4.0 X 7.5 X 0.75	10.1 x 19 x 1.9
INSTRUMENT - MINI TRAYS - 1600 SERIES			
TR-810-RO	Micro Tray - Base, Insert Tray, Lid & 2 Mats	6.0 X 10 X 1.5	15.21 x 25.4 x 3.8

NOTE: One year warranty from the date of purchase against failure due to manufacturing defects. Defective products will be replaced at no charge to the customer upon prepaid return of the product.