

## BATRIK ENDOSCOPY Reusable Channel Brushes DAYBRUSH™


*If a facility still prefers a Reusable brush, then the DayBrush™ is the perfect solution.*

BATRIK has developed a disposable channel brush that can be used for cleaning endoscopes for an entire day. You no longer need to worry if your present reusable brush is cleaning effectively. After processing ten (10) endoscopes<sup>(1)</sup>, simply discard the brush.

### The unique structure of the Reusable DayBrush™ enables the following

- soft tip ensures easy manipulation within the channel with no risk of damage
- the back-end is sealed and softened to prevent fluid from getting trapped inside the tube. The softening process also prevents damage if the cleaning brush is used in reverse
- can be steam sterilized
- the volume of bristles on the DayBrush™ is 20mm in length — almost double most other brushes. A large brush provides more surface contact during the cleaning process
- the BATRIK Reusable DayBrush™ has been tested to effectively clean up to 10 scopes<sup>(1)</sup>
- significant per use savings when compared to reusable or disposable channel cleaning brushes

<sup>(1)</sup> Discard after cleaning 10 scopes or sooner if any brush degradation is noticed. After ten applications the brush MAY have the appearance that it can still effectively clean, but we discourage prolonged use due to risk of cross contamination. Actual number of uses will vary based on handling and processing.

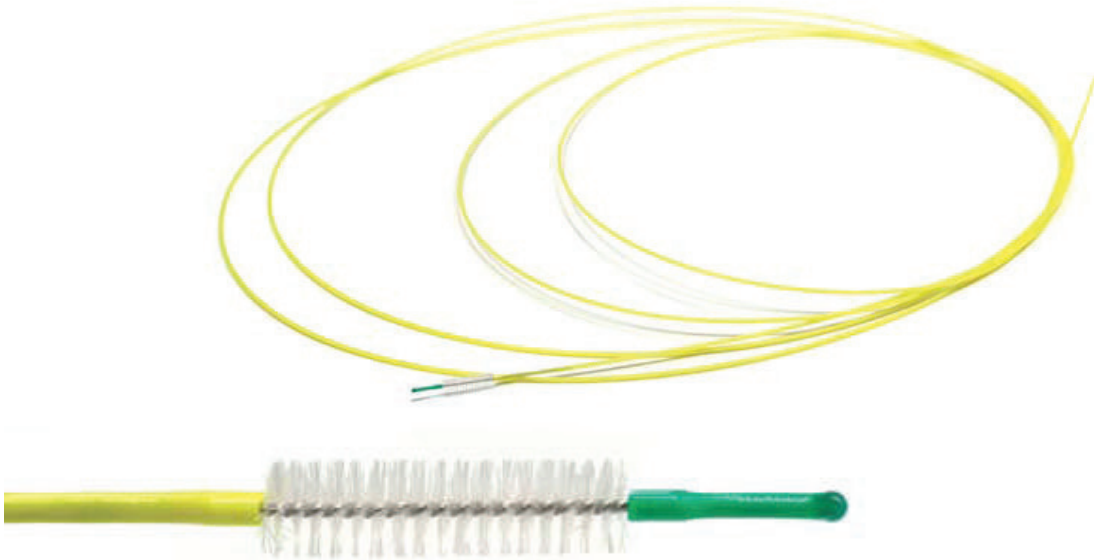

#### ADD THE SUFFIX...

F= receive a kit with our Flat Sponge (Packaging: Box of 40). Applicable on certain models only, call for details.

K= receive a kit with the single use Valve / Control Cleaning Brush 5/11 Double Brushes (item # 27-57787). Applicable on all models.

FK= includes both sponge and valve cleaning brush.



### 5mm Cylindrical DayBrush™

Item #	Length inches / cm	Brush Length inches / mm	Brush (Ø) Inches / mm	Sheath (Ø) Inches / mm	Brush Ends
30-18050-T	70" / 180cm	0.79" / 20mm	0.20" / 5mm	0.07" / 1.75mm	Single
30-23050-T	90" / 230cm	0.79" / 20mm	0.20" / 5mm	0.07" / 1.75mm	Single
30-23050-2T	90" / 230cm	0.79" / 20mm	0.20" / 5mm	0.07" / 1.75mm	Double

**Packaging:** 50 / Box