



**pfm**medical

# Proudly Serving Canadian Healthcare providers & patients

with a complete,  
universal drainage  
system, offering a  
standardized user-  
friendly & cost-  
effective one-stop  
solution.



**ASEPTI<sup>®</sup>**

Insertion Kits

Pleural & Peritoneal  
Drainage System

**ET-VA<sup>®</sup>**

Drainage Bottles

# SURGMED

GROUP 

Surgmed is dedicated to delivering innovative drainage solutions that support clinicians and improve patient outcomes across Canada.

Surgmed Group is a Canadian medical device distributor and manufacturer proudly serving the healthcare market since 1992 through its division, Advanced Surgi-Pharm. With decades of experience and a recognized presence in IPC (Intrapleural and Peritoneal Catheter) drainage nationwide, Surgmed partners with leading global technology providers to bring advanced, clinically proven solutions to hospitals, clinics, and specialty care teams.

Our experienced sales team — with dedicated representatives coast to coast in every province — provides highly on demand support for product selection, clinical training, and ongoing troubleshooting.

Through our strong partnerships with pfm, and hands-on IPC expertise, we ensure healthcare professionals have the tools they need to enhance outpatient care, improve patient comfort, and elevate quality of life.



## Why pfm?

pfmmedical is a leading medical device company specializing in interventional pleural and peritoneal drainage systems. Known for their innovation and clinically proven technologies, pfmmedical develops advanced solutions that help clinicians safely and effectively manage fluid buildup from conditions such as recurrent pleural effusion and malignant ascites. With a strong focus on improving patient outcomes and quality of life, pfmmedical's products are trusted by healthcare providers around the world and form a key part of Surgimed's portfolio of drainage solutions in Canada.

## ASEPT® Catheter and Insertion Kits

The ASEPT® Pleural and Peritoneal Drainage Systems are designed for long-term, intermittent drainage of symptomatic, recurrent fluid accumulation.



- |   |                                    |                        |                      |
|---|------------------------------------|------------------------|----------------------|
| 1 ASEPT® Catheter   | 6 16G x 1.0" (25.4mm) Needle       | 11 Foam Catheter Pad   | 16 10mL Syringe      |
| 2 ASEPT® Drainage Line  | 7 25G x 1.0" (25.4mm) Needle       | 12 Introducer Needle   | 17 J-Tip Guidewire   |
| 3 5-in-1 Drain Line Adapter                                       | 8 22G x 1.5" (38.1mm) Needle       | 13 Tearaway Introducer | 18 Tunneler          |
| 4 Chloraprep® 10.5mL Vial   | 9 Tegaderm® Self Adhesive Dressing | 14 8F (2.67mm) Dilator | 19 Forceps           |
| 5 Synthetic Suture, 2-0 Straight and Synthetic Suture, 3-0 Curved | 10 Gauze Pads, 4"x4"               | 15 12F (4.0mm) Dilator | 20 Safety Scalpel    |
|   |                                    |                        | 21 Fenestrated Drape |

The primary component of the ASEPT® Drainage system is the catheter, which has an aseptic, self-sealing valve attached that provides a secure and reliable connection. The valve can be connected to an ASEPT® vacuum drainage bottle, drainage bag, or a stand-alone drainage line to connect to other fluid collection equipment.

### Features:

- ▶ Secure connection
- ▶ Cap-less and Aseptic design
- ▶ Needle-free and swabable
- ▶ Easy to replace
- ▶ Self-sealing

## ASEPT<sup>®</sup> Catheter Portfolio

### Long-Term Intermittent Drainage Solution

The pfm ASEPT<sup>™</sup> Catheter System is designed to provide controlled, long-term, intermittent drainage for patients with pleural effusions and malignant ascites – supporting symptom relief in hospital or outpatient settings.

Here are a few key features of the PFM Catheter 15Fr

- Designed for long-term, intermittent drainage of pleural effusions (M7001-CAN) and malignant ascites (M7008-CAN).
- Aseptic, needle-free, swab-able valve helps maintain a secure and hygienic catheter connection.
- Carefully engineered to allow periodic access to the pleural or peritoneal cavity for symptom relief.
- Included in the kit with most items generally required for catheter implantation, streamlining the procedure.
- Works with compatible drainage lines, bottles, or bags for flexible clinical or outpatient use.

### Single Catheter Platform for Pleural & Peritoneal Use

Designed for procedural efficiency and flexible patient management.

#### 15Fr Silicone Catheter

Smooth, flexible design for patient comfort and reliable, long-term placement.

#### M7001-CAN - 66 cm

Designed for pleural effusions.

#### M7008-CAN - 71 cm

Designed for malignant ascites



#### Aseptic Needle-Free Swabable Valve

Helps maintain a secure, hygienic connection and reduces infection risk.

#### Engineered for Intermittent Drainage

Allows periodic access to the pleural or peritoneal cavity for symptom relief.

#### Complete Implantation Kit

Includes most items generally required for safe and efficient catheter placement.

**ASEPT® Roadmap**

Integrated Pathway for Long-Term Intermittent Drainage

**1**

**Patient Consultation & Catheter Placement**

Patients are scheduled and seen by their surgeon for ASEPT® catheter insertion either Pleural (M7001-CAN) or Peritoneal (M7008-CAN). Patients equipped with an IPC catheter require ongoing drainage as directed by their healthcare provider.



**2**

**pfm Universal Valve Replacement**


Many patients have IPC catheters with brand-specific valves that are only compatible with certain bottles. This creates challenges during urgent or emergency drainage situations. pfm offers a universal Luer-Lok valve that can be replaced bedside, allowing compatibility with any Luer-Lok connection (including wall suction), and most importantly, enabling use with pfm drainage bottles for immediate relief.



**3**

**Clinical Drainage - Jet-Vac**

Following insertion, patients may require periodic drainage in clinic. pfm Jet-Vac vacuum bottles are available in 600 mL, 1000 mL, and 1500 mL volumes. Options include attached drain lines, valve-only configurations, and sampling ports for patients requiring additional testing.

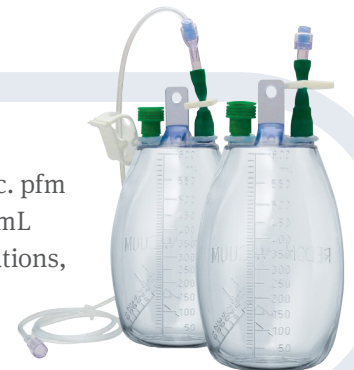

**4**

**ASEPT® Patient Homecare**

For patients draining at home, pfm offers ready-to-use, patient-friendly drainage kits. These include personal-use bottles and all necessary dressing materials in a convenient carrying case.

This solution is ideal for both individual patient use and home healthcare teams conducting home visits.



Designed for long-term, intermittent drainage of symptomatic fluid accumulation.

# ET-VA Drainage Bottles

## Breakdown

- ▶ Available in three sizes:
  - ▶ 600 mL
  - ▶ 1,000 mL
  - ▶ 1,500 mL
- ▶ Available in two configurations:
  - ▶ Drain line featuring sampling port for easy access to contents of bottle via luer-lock syringe:
    - Connection is made at catheter level
    - Drain line cannot be removed or changed
  - ▶ Luer-lock valve:
    - Connection is made at bottle level
    - Fluid samples can also be drawn with luer-lock syringe connected directly to valve
    - Internal Sterility



## ASEPT®

### Breakdown

- ▶ Available in two sizes:
  - ▶ 600 mL
  - ▶ 1,000 mL
- ▶ Only available with drain line can be used to regulate flow, should patient feel any discomfort during drainage procedure.

## Differences between ASEPT<sup>®</sup> and Jet-Vac Drainage Systems





Jet-Vac Centesis Catheters, Jet-Vac Procedure Trays and Jet-Vac Bottles provide a complete drainage system for physicians, nurses, and hospitals. The Jet-Vac drainage bottle is a disposable vacuum device used to aspirate, remove, or sample body fluids.

## Jet-Vac Centesis Catheters

**Large circular side holes designed to increase fluid drainage**

**Tapered catheter tip**

- Provides a smooth transition
- Helps minimize patient discomfort

**HDPE material**

- Reduces catheter collapsing
- Promotes a smooth catheter insertion

**Clear needle hub with catheter lock**


- Increases fluid visibility
- Ergonomic design provides a secure grip
- Prevents needle/catheter separation

**Plastic vent plug**

- Allows for air to exit
- Provides added safety for physician

**Echogenic needle**

- Ensures needle visibility during ultrasound placement

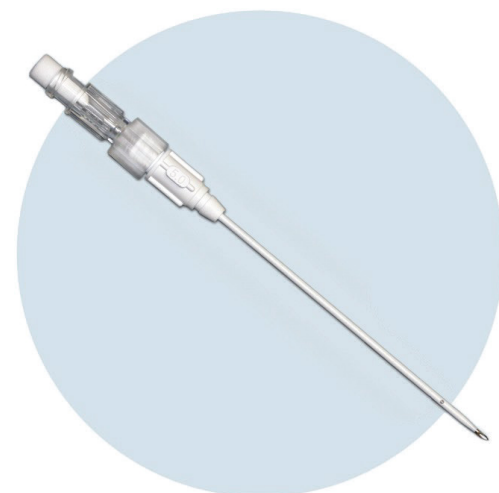

## Ordering Information

### Jet-Vac Centesis Catheters

Item No.	Catheter Size	Length	Needle Gauge	Box Qty
M7511	4F	7 cm	20G	5/Box
M7512	4F	10 cm	20G	5/Box
M7513	5F	7 cm	19G	5/Box
M7514	5F	10 cm	19G	5/Box
M7515	5F	15cm	19G	5/Box

**ASEPT® Accessories**

pfm medical, inc. is proud to offer a full line of accessories ensuring that the physician or patient will have everything he or she needs to complete the drainage procedure. Our accessory portfolio contains the following items:

ASEPT® Replacement Valve



ASEPT® Drainage Line Set



ASEPT® Procedure Pack



ASEPT® Permanent Drainage Catheter



Universal Connection to Drainage Bottle



**Ordering Information**

**ASEPT® Drainage Bottles and Bags - Patient Home Care**

Item No.	Description	Box Qty
M7002	ASEPT® 600mL Drainage Kit (Includes vacuum bottle and Procedure Pack)	10 each
M7000	ASEPT® 1000mL Drainage Kit (Includes vacuum bottle and Procedure Pack)	10 each
M7006	ASEPT® 600mL Select Drainage Kit (Includes low vacuum bottle and Procedure Pack)	10 each
M7007	ASEPT® 1000mL Select Drainage Kit (Includes low vacuum bottle and Procedure Pack)	10 each
M7060	ASEPT® 2000mL Drainage Starter Kit (Includes gravity bag with drainage line - for Peritoneal drainage only)	10 each
M7058	ASEPT® 2000mL Drainage Bag (Includes gravity bag with drainage line - for Peritoneal drainage only)	10 each
M7080	ASEPT® 600mL Drainage Starter Kit (Includes vacuum bottle and Procedure Pack in patient box)	5 each
M7070	ASEPT® 1000mL Bottle (Includes vacuum bottle - does not include Procedure Pack)	10 each
M7072	ASEPT® 600mL Bottle (Includes vacuum bottle - does not include Procedure Pack)	10 each
M7081	ASEPT® 1000mL Drainage Starter Kit (Includes vacuum bottle and Procedure Pack in patient box)	5 each

## Ordering Information

### ASEPT® Insertion Kits

Item No.	Description	Box Qty
M7001-CAN	ASEPT® Pleural Drainage System (Includes ASEPT® catheter and insertion kit)	1 each
M7008-CAN	ASEPT® Peritoneal Drainage System (Includes ASEPT® catheter and insertion kit)	1 each
M7009	ASEPT® Supplemental Insertion Kit (Includes 7F dilators & 16F valved peelable introducers)	10 each

### Jet-Vac Bottles - Clinical Use

Item No.	Item	Description	Box Qty
M7090	Jet-Vac 1000	1,000 mL Evacuated Drainage Bottle in Valve-Only Configuration	10/Case
M7092	Jet-Vac 600	600 mL Evacuated Drainage Bottle in Valve-Only Configuration	10/Case
M7093	Jet-Vac 1500	1,500 mL Evacuated Drainage Bottle in Valve-Only Configuration	10/Case
M7094	Jet-Vac 1000D	1,000 mL Evacuated Drainage Bottle with Drain Line	10/Case
M7095	Jet-Vac 1500D	1,500 mL Evacuated Drainage Bottle with Drain Line	10/Case
M7096	Jet-Vac 600D	600 mL Evacuated Drainage Bottle with Drain Line	10/Case

### Accessories

Item No.	Description	Box Qty
M7004	Catheter Replacement Valve	5/Box
M7005	Drainage Line Set (Includes Drainage Line with End Caps, and 5-in-1 Adapter)	10/Case
M7040	Procedure Pack (Traditional Dressing Change)	10/Case
M7097	Jet-Vac Drainage Line Extension (two female luer connectors)	10/Case
M7098	Jet-Vac Large Bore Drainage Line	5/Box
M7099	Jet-Vac Drainage Line	10/Case
M7555	2000cc Drainage Bag	10/Case

**ADVANCED**  
SURGI-PHARM INC.

Scan the QR Code to learn  
more about our products



Pleural Peritoneal  
Drainage System



Jet-Vac Bottles