

introducers

For PICC catheters

Stand-alone introducers for neonatal PICC catheters: splitting needle, microflash® peelable cannula and IV cannula.



IV cannula

IV cannula

Product Code	Size	For catheter	Units/Case
CMS-207	24 G	1 Fr	10



Splitting needle

Splitting needle

Product Code	Size	For catheter	Units/Case
1261.20N	24 G	1 Fr	20
1262.30N	20G	2 Fr	20




microflash™

Peelable Introducer/Peelable Safety Introducer*

Peelable Introducer for the insertion of PICC & Midline catheters, with notched needle and clear cannula sheath for early detection of blood flashback

Product Code	Size	For catheter	Units/Case
7371.19	20 G	1 Fr / 2 Fr	10

- A** Top hole needle shaft permits blood flashback between needle & cannula for immediate indication of vessel access
- B** Short, lightweight safety housing for balance and ease of handling
- C** Curved hub for secure ergonomic fingertip placement
- D** Tri-bevel needle tip for ease of skin penetration
- E** Smooth needle to cannula transition reduces drag through the tissue
- F** Transparent cannula permits immediate blood visualization upon vessel penetration
- G** Funnel cannula opening to facilitate easy catheter entry
- H** Radiopaque stripe on cannula for patient safety

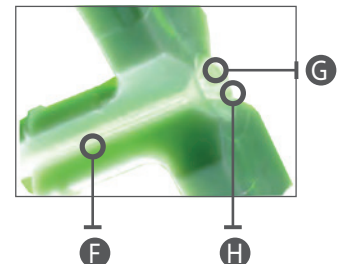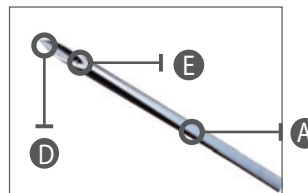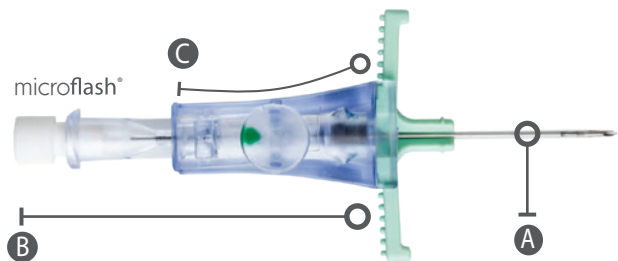

*The microflash® Introducer is intended for use as a split cannula introducer used to puncture peripheral veins.